

Get Started with Q-interactive®

Q-interactive Start-Up Guide

Here are instructions for installing the Q-interactive Assess app, establishing the proper settings on both iPad® devices, finalising your login credentials, and creating a practise assessment session. We recommend practising the tests a few times before assessing real clients. Sign-up to one of our monthly webinars to get more information about the setup and use of Q-interactive for new users.



This is a good opportunity to ask any questions about Q-interactive that you may have. Visit: pearsonclinical.co.uk/webinars to sign-up.

SECTION 1: Download the App


The Q-interactive Assess app must be downloaded onto both iPads in order to use the system.

Step 1



Open the App Store, found on the home screen of the iPad.

Step 2

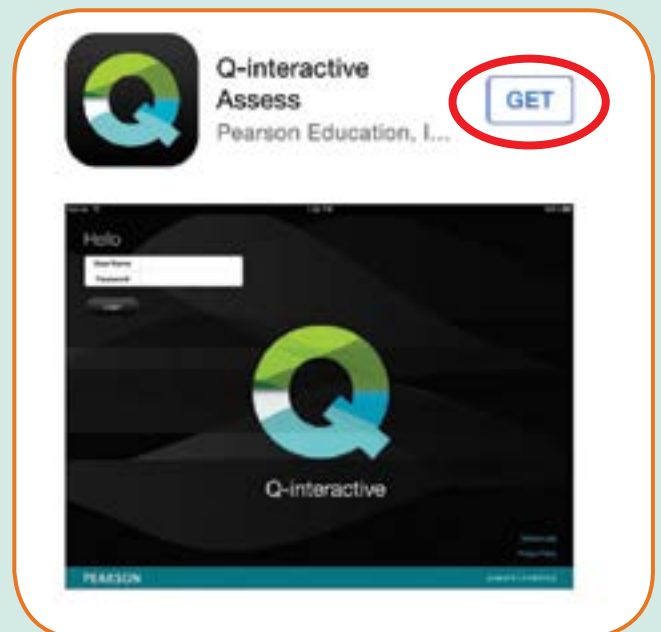

Search for Q-interactive and tap on "GET."

SECTION 2: Configure the iPads

With Assess downloaded, let's configure the settings on both devices. This is an essential step for a few reasons:



- Ensure iPads pair correctly in order to capture data from each other.
- Stimuli display in a standard manner.
- Assist in the smooth operation of Q-interactive.
- Eliminate distractions, which may occur during the assessment.

Step 1



Touch the Settings app located on the home screen.

Step 2: Adjust your iPads according to the following table.

Practitioner's iPad	Client's iPad	Setting	Menu	Choice
•	•	Wi-Fi	On	Choose network
•	•	Bluetooth	On	
	•	Display & Brightness	Brightness	(slide bar to far right)
•	•	Display & Brightness	Auto-Lock	15 Minutes
	•	General	Accessibility > Magnifier > Auto-Brightness	
•		Assess	Application Mode	Practitioner
	•	Assess	Application Mode	Client

SECTION 3: Creating an Assessment Session

Your Welcome email gets this started by first establishing your login credentials.


Step 1



The "Click here" link will open the screen to establish your Q-interactive password.

After submitting the password and accepting the Terms and Conditions of Q-interactive, you are now ready to create a practise assessment session.

Step 2

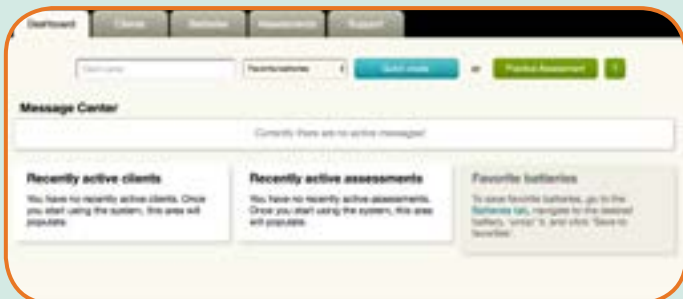

From the Dashboard tab, click the Practise Assessment button to make the Practise Client profile and Save.

Step 3



Select battery by clicking the button to the left, then click the Next button.

Step 4



Set the date and time for administration and click the Finish button.


Step 5



Click the Send to iPad button and you'll be all set to log into the iPads!

Login to both iPads with the same login credentials established earlier to begin testing.

Step 6

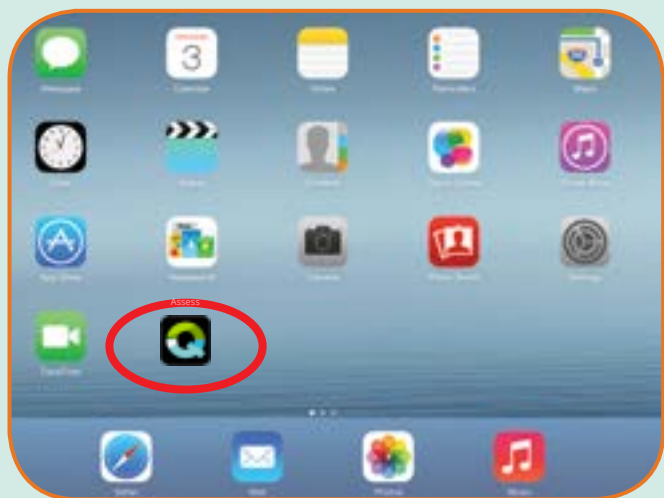

Step 7



Enter your User Name and Password. Allow 10-20 minutes for your complete library of tests to download.

Your iPads will connect when you start your assessment and both iPads have been configured correctly.

Step 8



Tap the Battery tab, then the Start button to begin the assessment. Practise assessment sessions have all the functionality of real assessment sessions, except all scaled/standard scores default to zero.

Congratulations! You are now ready to use Q-interactive!

Please reference the Q-interactive User Guides and training videos found on the Support tab within Q-interactive Central (qiactive.com) for detailed instructions. qiactive.com for detailed instructions.

0345 630 8888* | pearsonclinical.co.uk/qisupport

Copyright © 2018 Pearson Education, Inc. or its affiliates. All rights reserved. Pearson and Q-interactive are trademarks, in the U.S. and/or other countries, of Pearson plc. CLINA15779 - 15739 SR 9/18

*Calls cost 3p per minute plus your phone company's access charge.

